
DentalDeal


License No.: 20180203 in Hangzhou,Zhejiang
Standard No.: QZJXB04-2009

ALGINATE MIXER INSTRUCTION

HL-YMC Ⅰ HL-YMC Ⅱ

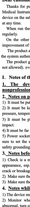

HL-YMC Ⅲ HL-YMC Ⅳ

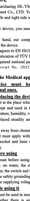

Hangzhou ZhongRun Medical Instrument Co., LTD.

Dear customer:

Thanks for purchasing this HL-YMC Alginate mixer produced by Hangzhou ZhongRun Medical Instrument Co., LTD. You must read the instruction carefully before using the mixer so that we will get safe and stable mixer to be used for a long time and to be at no loss and at any time.

When use the device, you must observe the rules strictly and preserve correctly and safely.

On the other hand, our company has the right to reserve an extra notice of the improvement of our products.

The product agrees to EN ISO13485: 2003+AC 2009 standard. which is approved by the centre accreditation of YUE Metrological Institute.

The product passed national patent: patent No: ZL 2007 2 0106233.0 (constructing patent) and patent No: ZL 2007 3 0106234.5 (appearance patent).

1. Notes of the Medical appliances

1) The device must be used to specific persons instead of nonprofessional users.

2) Before operating the device:

1. The mixing tube should be well connected with vessel or any other liquid.
2. It must be kept and used in the place where there is no adverse air caused by air pressure, temperature, humidity, ventilation, dust, sunshine, and salt, ion, etc.
3. It must be placed steadily and stably without any tilt, quake or the possibility of tilt.
4. It must be far away from chemicals and chemical gas.
5. Power socket must apply with the required voltage and current of the device. Make sure that the socket and lines connected with the device is reliable grounded to avoid the safety grounding.

3) During using the device:

1. Check up the work before using. Make sure there is no abnormal phenomenon on the appearance, eg: on worm, the electrical lines of the socket is integrated without any crack or breakage, the vertical and adjacent is correct.
2. Make sure the safety grounding correctly and stably.
3. Make sure the adjusting voltage is the same as power with the requirement.

4) Notice while using it

1. The device must be used continuously, eg it used for mixing impression material.
2. Monitor whether there is any abnormal phenomenon, eg: noise, smoke. Once it is abnormal, turn off the switch and pull out the socket.
3. When the mixing tube is blocked, you can take the internal force is too big to overloading the motor, the machine will stop working automatically within 5 seconds, this is a kind of protection. Please turn off the power button, or stop on the foot control switch to pause to start working it.
4) The patient mouth's context the device.

5) The patient mouth's glasses to avoid spatters.

5. Notes after using

1) Checks must be made according to the requirement strictly and turn off the switch.
2) After turning off the switch, the socket must be held down instead of pulling down.
3) Clean the mixing tube and the connecting tube.
4) Clean the device with a soft cotton cloth, and wipe it with a dry soft cotton cloth.
5. It must be kept and used in the place where there is no adverse air caused by air pressure, temperature, humidity, ventilation, dust, sunshine, salt, ion, etc.
6. It must be placed steadily and stably without any tilt, quake or the possibility of impact.
7. It must be far away from chemicals and chemical gas.
8. The device must be cleaned and kept dust-free. The socket and its accessories must be kept dust-free and dry.

6. You power if dismantle, or repair the device by yourself when not provision, must context distributors, and repair be specific persons.

7. You mustn't remake the device, or you have to be responsible for the consequences and lose the warranty at last.

8. Regular check and maintenance

1) Under normal circumstances, it is suggested that check once before using and make maintenance.

2) If it is in idle for a while, the check on the safety and function is needed before operation to avoid the accidents as far as possible.

9. Caution

1. Pay attention to spatula and mixing bowl when operation, avoiding water and unnecessary material spill out to hurt somebody or something.

2. Because there may occur induction electric or leakage, the device must be far away from the patient's eyes.

3.The device should be cleaned and placed in dry and ventilated room without corrosive gas, and avoid the direct sunlight.

4.To prevent the water and gum from entering the motor, the lid must be kept tight and used moist phases.

5. The instruction explains the safe use of the device specifically, please be sure to put it near the device and read it carefully and frequently.

6.Your company is not responsible for the following problems:

1. Mechanical damage and personal accident caused because the staff out of company's training operate the device rules.

2. Mechanical damage and personal accident caused due to using the parts out of our recognition.

3. Mechanical damage and personal accident caused because of not obeying the above rules promptly.

4. Mechanical damage and personal accident caused due to violating the required voltage, grounding and placed environment.

5. Mechanical damage and personal accident caused by natural disasters such as earthquake, fire, etc.

6. Spatula and mixing bowl are easy to be damaged, the fee on renewing and changing them should be dealt by yourself.

Marks on the introduction:

Caution

Mechanical damage and personal accident will be caused if you violate it.

Attention

Mechanical damage and personal accident will be caused if you violate it.

III. Type and Specifications

1. **Product Name**: Alginate mixer

2. **Specifications**: HL-YMC1
HL-YMC II
HL-YMC III
HL-YMC IV

3. **Power supply**: AC 110V, 60 Hz ①
AC 120V, 60 Hz ②
AC 220V, 50 Hz ③
AC 230V, 50 Hz ④
AC 240V, 50 Hz ⑤

4. **Motor power**: 45 W

5. **Comprehensive performance comparison**

TYPE	SPEED	Timing Range /+/-	WEIGH(Kg)	SIZE(mm)
HL-YMC I	3000rpm/2800rpm (high-speed/-20%rpm)	Regular alarm/-20s, no timing	2300	260*220*190
HL-YMC II	3000rpm/2800rpm (high-speed/-20%rpm)	Regular alarm/-20s, no timing	2360	260*220*190
HL-YMC III	3000rpm	Regular alarm/-20s, no timing	2360	260*220*190
HL-YMC IV	0-3000rpm	Regular alarm/-20s, no timing	2360	260*220*190

6. **Mixing Time**: according to the specification of the impression material manufacturer.

7. **Classification of the devices**: Ⅱ-B: Electrical devices. common devices. non -cp# or non-APG devices, use ms environment.

8. **Service life**: 5 years.

IV. Names of parts and their specifications:

HL-YMC I



1. Shell 2. Mixing bowl 3. Power supply indicator light 4. Spatula
5. Timing indicator light 6. Power switch 7. Base 8. Power cord 9. Nameplate
10. Start/Start button (low speed/high speed changer) 11.Stop button

HL-YMC II



1. Shell 2. Mixing bowl 3. Power supply Indicator light 4. Spatula
5. Timing indicator light 6. Power switch 7. Base 8. Power cord 9. Nameplate
10. Start/Start button (low speed/high speed changer) 11.Stop button 12.Foot switch
13. Foot switch connector

Operating Instruction: HL-YMC I

HL-YMC II

After 20S working, the warning light will flash and sound a beep

START SET:
First press is low speed, second press is high speed



① Power switch: a black boot-shaped switch behind the machine.

② Low speed: Press the "Start" button, the mixing bowl starts rotate, the mixing bowl turn in the low speed. At the same time, the mixing bowl can unloaded.

③ "Start/Set" button: the "Start/Set" button in a left side of the front panel. The first press starts to the low speed; it will become high speed after 20S.


④ "Stop" button. The "Stop" button is at right side of the front panel, press "Stop" button and the working immediately.

HL-YMC III




1. Shell 2.Mixing bowl 3. Power supply Indicator light 4. Foot switch
5. Spatula 6. Warning indicator light 7. Power switch
8. Foot switch connector 9. Fuse 10. Power cord 11. Nameplate

HL-YMC IV

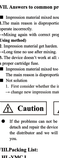

1. The names of parts are the same as HL-YMC III.
2.Each function display icon

Operating Instruction: Indicator diagram of the foot switch connect with the mainframe.

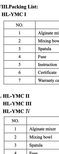

V. Usage:

1. HL-YMC I

1. Put the cleansing and disinfected mixing bowl into the tray, turn clockwise to lock it.
2. Let power plug connect to the electric socket with grounding, turn on the power switch.
3. Put the water and impression material into the mixing bowl according to the proportion from the manual of impression material. Then mix to blob state.
4. Put the spatula into mixture, press Start/Set button(S) to start at low speed (the liquid will not splash out), after 20 seconds working, the warning sign flashes and sounds, the spatula stirring the wall of mixing bowl, it's better for mixing, after 20S the alarm light will flash, then to notice user to use the mixture as soon as possible.
5. Press Start/Set button (S) again to stop working. Use spatula along mixing bowl to remove the mixture, just complete.
7. Anticlockwise rotate the mixing bowl, then take it out, and clean it.

2. HL-YMC II

1. Put the cleansing and disinfected mixing bowl into the tray, turn clockwise to lock it.
2. Let power plug connect to the electric socket with grounding, turn on the power switch.
3. Put the water and impression material into the mixing bowl according to the proportion from the manual of impression material. Then mix to blob state.
4. Put the spatula into mixture, press Start/Set button(F) to start at low speed (the liquid will not splash out), after 20 seconds working, the warning sign flashes and sounds, the spatula stirring the wall of mixing bowl, it's better for mixing, after 20S the alarm light will flash, then to notice user to use the mixture as soon as possible. Then it turns to high speed.
5. Press Stop button (T) to stop working. Use spatula along mixing bowl to remove the mixture, just complete.
6. Anticlockwise rotate the mixing bowl, then take it out, and clean it.

3. HL-YMC III

1. Put the cleansing and disinfected mixing bowl into the tray, turn clockwise to lock it.
2. Let power plug connect to the electric socket with grounding, turn on the power switch.
3. Put the water and impression material into the mixing bowl according to the proportion from the manual of impression material. Then mix to blob state.
4. Put the spatula into mixture, step foot control lightly to start at low speed (the liquid will not splash out), after 20 seconds working, the warning sign flashes and sounds, the spatula stirring the wall of mixing bowl, it's better for mixing, after 20S the alarm light will flash, then to notice user to use the mixture as soon as possible. Then it turns to high speed.
6. Step foot control, it stroke to automobile accelerator.
7. Anticlockwise rotate the mixing bowl, then take it out, and clean it.

4. HL-YMC IV

1. Put the cleansing and disinfected mixing bowl into the tray, turn clockwise to lock it.
2. Let power plug connect to the electric socket with grounding, turn on the power switch.
3. Put the water and impression material into the mixing bowl according to the proportion from the manual of impression material. Then mix to blob state.
4. Put the spatula into mixture, step foot control lightly to start at low speed (the liquid will not splash out), after 20 seconds working, the warning sign flashes and sounds, the spatula stirring the wall of mixing bowl, it's better for mixing, after 20S the alarm light will flash, then to notice user to use the mixture as soon as possible. Then it turns to high speed.
5. Step foot control, it stroke to automobile accelerator.
6. Adopt the large LCD display, which can bring vivid image, let you be full of wit and humor in operation.
9. Leave foot control, the machine will stop. Use spatula along stirring bowl to remove the mixture.
10.Anticlockwise rotate the mixing bowl, then take it out, and clean it.

Attention

The mixing time should not be too long, (generally within 20S), the effective using time according to the manual of impression material manufacturer. After well mixing, please take out mixture promptly to avoid dry condensate.
3. After using, must turn off the switch and confirm the indicator light is off.

The notice during using:

Attention

A. Long time mixing may result in bad performance of impression material.
B. After using, must turn off the switch.
C. Please unplug the power cord if you do not use in a long time.

warning

A. The equipment is not water-proof.
B. Do not put the mixture into motor, please don't use it for other materials, especially toxic liquid and chemicals.
C. Please don't touch the motor parts when using to avoid the harm of the human body. Especially the metal part between the bowl and body.

VI. Daily inspection and maintenance

A.Turn off the power switch.

warning

1) The inspection and maintenance must be done after complete cooling.
2) The plug must be pulled down.
3. Please clean the bowl and spatula each time after using.
4) Wipe away the dust on surface.

Attention

Chemical cleaner mustn't be used for cleaning the device.

VII.Answers to common problems

■ Impression material mixed nonuniform or pat nonuniform.

1.The main reason is disproportion of water and impression material, mix too short or too long time.

2.Mixing speed is incorrect proportion, correct operation.(Details please check 5VII. **Using method**)

3. Impression material get harden.

1)Long time to use after mixing.

2)The mixing tube is not at all clock whether the fuse is broken, if it is, please change a same amperage fuse.

■ Impression material mixed too sparse or too thick.

The main reason is disproportion of water and impression material.

■ Not solution

1. First consider whether the impression material is deteriorated.
2. change new impression material and operate it again

Caution **Attention**

● If the problems can not be solved according to the above suggestions, don't operate and repair the device by yourself. It is a must to keep in connection with the distributor and we will supply our specialized staff to solve them for you.

VIII.Packing List:

HL-YMC I

NO.	Name	Quantity(piece)
1	Alginate mixer	1
2	Mixing bowl	1
3	Spatula	1
4	Fuse	1
5	Instruction	1
6	Certificate	1
7	Warranty card	1
8	Foot switch	1

2. HL-YMC II

HL-YMC III

HL-YMC IV

NO.	Name	Quantity(piece)
1	Alginate mixer	1
2	Mixing bowl	1
3	Spatula	1
4	Fuse	1
5	Instruction	1
6	Certificate	1
7	Warranty card	1
8	Foot switch	1